

Specifications

- Battery: Lithium Ion, DC 3.6-4.2V
- Power adapter: 90-250V AC, 50-60 Hz
- Class: Class II (ICE601-1)
- Charging base: 9V DC
- Power consumption: less than 10W
- Wavelength range: 430-490 nm
- Light Intensity: 1,250 mW / cm²
- Dimensions: 7"L x 1"W x 1"H
- Weight: 160 g
- Warranty: 2 years



Benefits

- Ergonomic, light weight design
- High intensity output up to 1 hour before recharging
- Cordless for ultimate convenience
- Cures most materials in 5-10 seconds
- Curing time: 20 seconds per cycle with beep every 10 seconds
- Additional battery included

Charging Stand Light Indicators

LED color	Light Strength
Green	Normal
Orange	Light strength insufficient

Curing Status Light Indicators

LED color	Status Indication
Green	Normal
Red	Battery Voltage Low
Red Flashing	LED curing light handpiece is over temperature

Curing Mode Light Indicators

LED color	Working Mode
Green	Fast curing: Time setting 20 seconds (every 10 second, it sounds once)
Red	Pulse curing: interrupt radiating (every 10 seconds, it sounds once)
Orange	Step curing: soft start for 10 second from 0 to full power, then full power during the later 10 second






Charging Stand



Curing Light



Model LED 55N

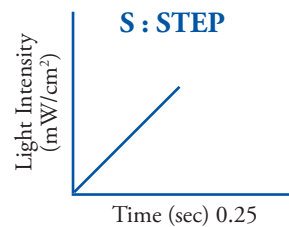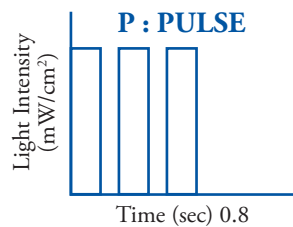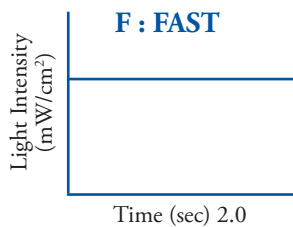

Shipping Dimensions: 10"L x 9"W x 5"H / 3 lbs.
(10/cs.) 25"L x 11"W x 17"H / 25 lbs.