
Unit head plumbing schematic

TPC

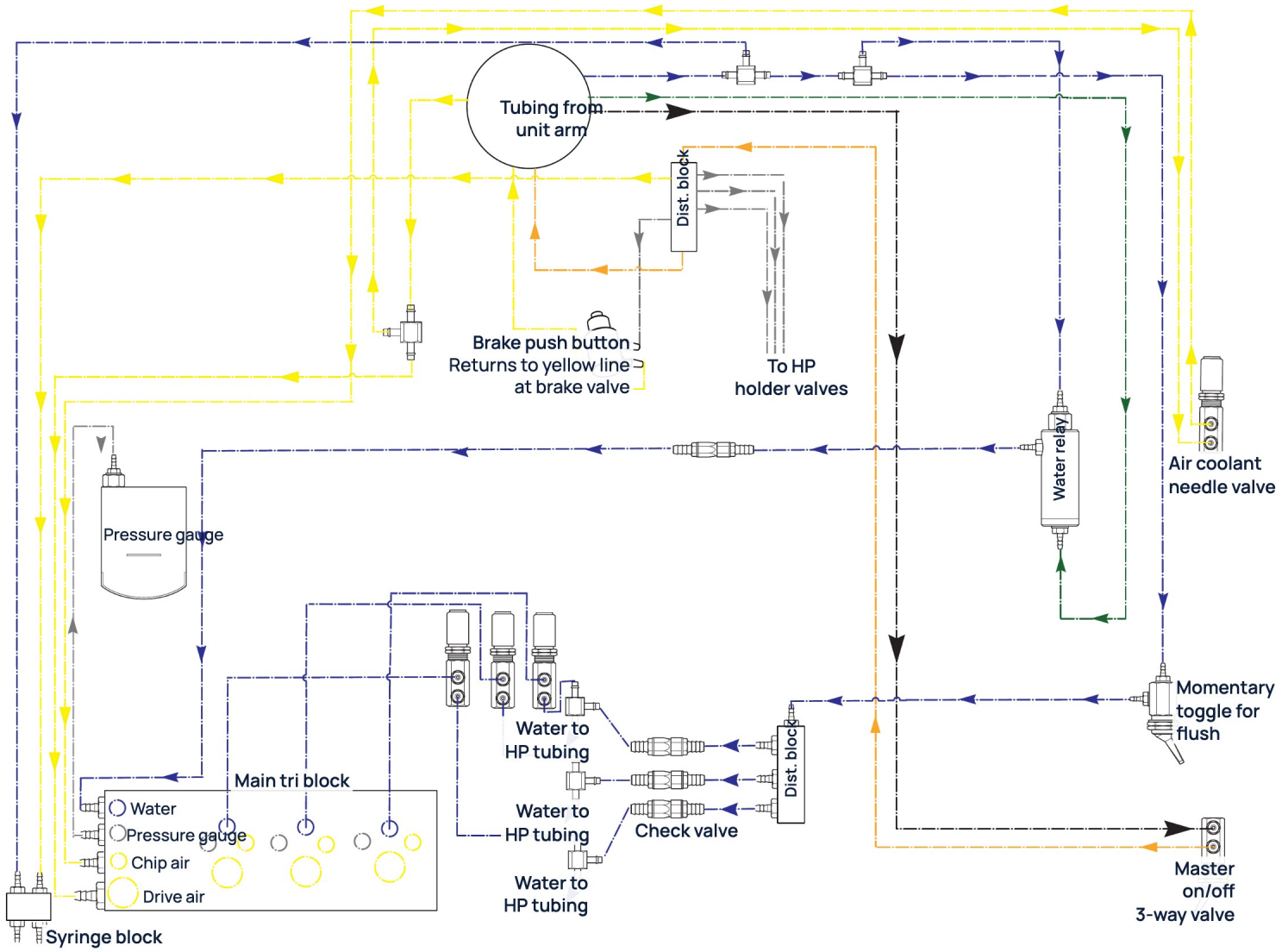

LINE KEY

- **Water**
- **Air** · chip and drive air
- **Supply air** to master on/off valve
- **Return air** from master on/off switch
- **Signal air** for the wet/dry toggle
- **Pilot air** to HP holder valves and pressure gauge

NOTES

Arrows show direction of flow.

Tri block circles match line colors: blue water, gray pressure gauge, yellow chip and drive air.


Drawing not to scale.

See sheet 2 for side box and junction box plumbing.

DRAWING	MODELS	REVISED	SHEET	TPC DENTAL
Mirage 2.0 unit head plumbing	2015 / 2000 (2.0)	October 2025	1 of 2	(800) 560-8222 · tpcdental.com

Side box and junction box plumbing schematic

TPC

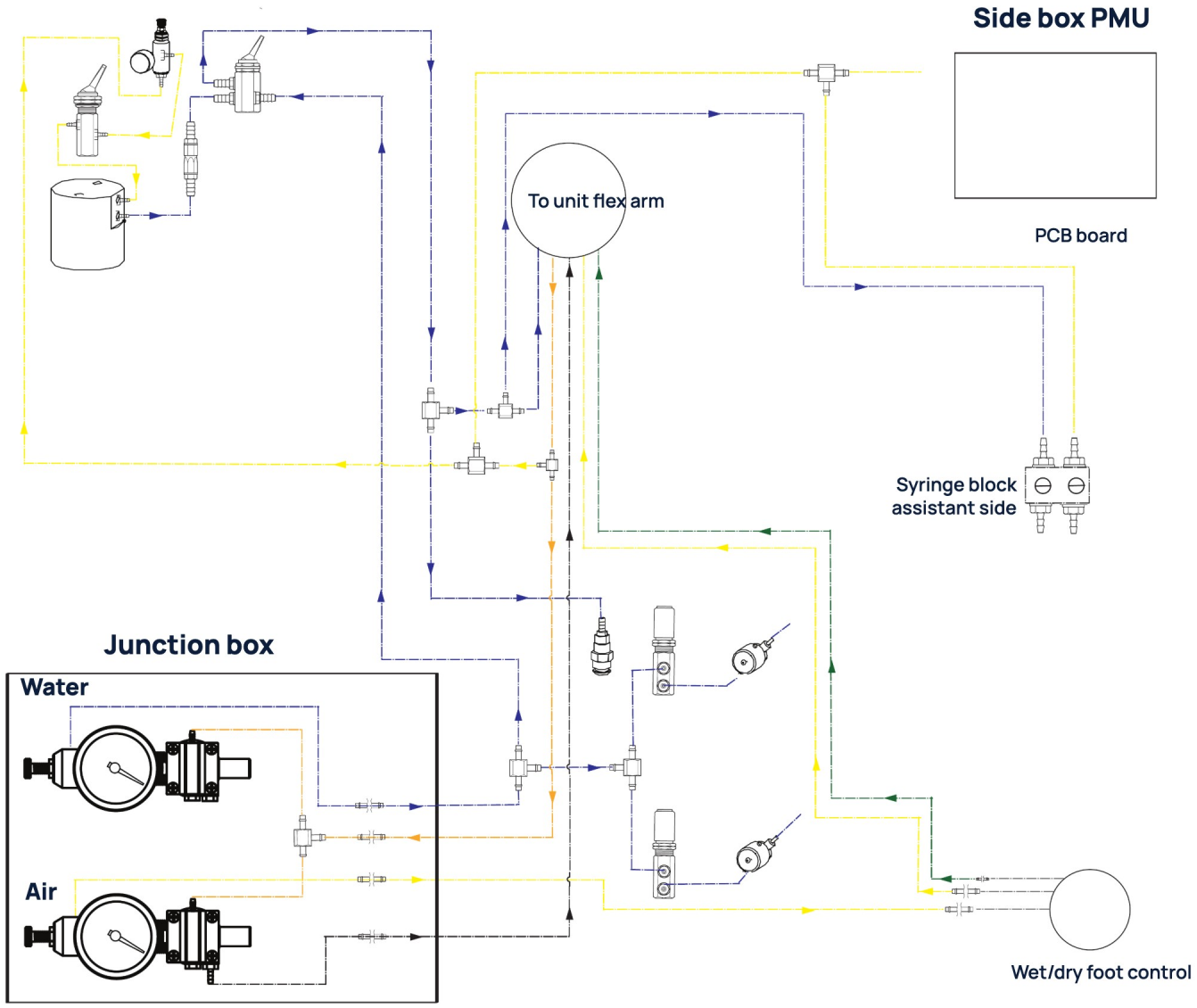

LINE KEY

- **Water** from junction box water regulator
- **Air** from junction box air regulator
- **Supply air** to master on/off valve
- **Return air** from master on/off switch
- **Signal air** for the wet/dry toggle

NOTES

Arrows show direction of flow.

The junction box holds the water and air regulators with gauges. All supply lines run from here to the side box, flex arm and foot control.

Drawing not to scale.

See sheet 1 for unit head plumbing.

DRAWING	MODELS	REVISED	SHEET	TPC DENTAL
Mirage 2.0 side box plumbing	2015 / 2000 (2.0)	October 2025	2 of 2	(800) 560-8222 · tpcdental.com