
Mirage 2.0 Floor Template

TPC

PLAN AND SIDE SECTION

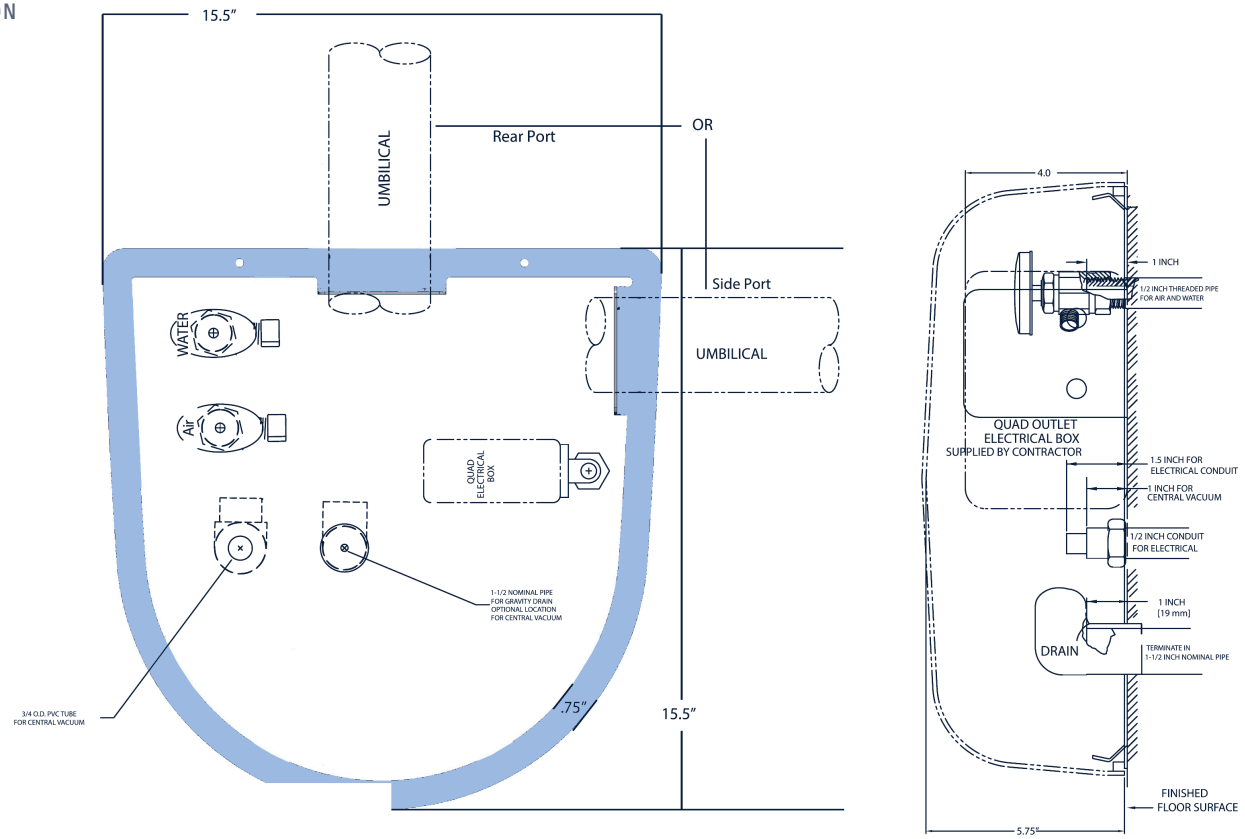

UTILITY CENTER, W AND D	WALL	SECTION DEPTH	SECTION BASE	AIR / WATER STUB
15.5″ 394 mm	**0.75″** 19 mm	**4.0″** 102 mm	**5.75″** 146 mm	**1″** 25 mm

Air and water ½″ threaded pipe **Electrical** 1½″ conduit, ½″ conduit to box **Central vacuum** 1″ or ¾″ O.D. PVC tube **Drain** 1½″ nominal pipe **Umbilical** rear or side port

1. Local regulations require licensed plumbers and electricians to install utilities.
2. All plumbing must conform to prevailing local codes.
3. Central vacuum plumbing to the utility center is specified by the vacuum supplier and terminates in the utility center with a ¾″ O.D. tube, perpendicular to the floor like the drain.
4. The gravity drain is a 1½″ nominal pipe protruding from the finished floor. Place a trap in line per local codes (supplied by contractor).